

Valiant 2 rehab XXL

Robust and reliable treadmill for rehabilitation



Highlights

High standards

Lode is a socially and environmentally responsible company. All Lode products are RoHS/WEE compliant and Lode is ISO 9001:2015, and ISO 13485:2016 certified. All medical products comply to MDD 93/42/EEC, incl. IEC 60601-1.

Various handrails available

Various handrails are optional available, making the treadmill suitable for your specific stress test setting.

Smooth acceleration

All treadmills have a smooth acceleration from 0 km/h to the start speed.

Extremely accurate

The treadmill is extremely accurate in its speed and angle settings

Patient friendly treadmill

The treadmill offers ultimate comfort for the patient:

- low step-up height
- faultless operation
- smooth acceleration
- reliable and reproducible test results
- low noise



Valiant 2 rehab XXL

Robust and reliable treadmill for rehabilitation



Modern designed treadmill specifically designed for rehabilitation. The low step-up height makes the treadmill perfectly suitable for all test subjects in the field of rehabilitation. The emergency stop with magnetic lanyard on the standard front handrail provides additional safety for the user. The Valiant 2 rehab offers a smooth acceleration from 0 km/h and is continuously adjustable in a range of 0,1 – 20 km/h. The network module provides control within a Lode software network (LEM, LRM or LCRM). Thanks to the unique low design of the cover plate it is possible to place a mirror or camera in front of the treadmill which makes it easier to monitor the gait of the test subject. The Valiant 2 rehab XXL has a standard running surface of 250 x 120 cm and 0-15% elevation. The treadmill is supplied with a 7" control unit with touchscreen for manual control and a bottle holder.

A USB A-B cable only for service purposes will be standard delivered with the product.

To connect LEM, LRM or LCRM you need a special interface cable that can be ordered under part number 930930.

Features



Low noise

Material choices, refined components and accurate manufacturing techniques lead to low noise.



Robust design

The product is designed to withstand continuous heavy use by subjects in most weights and sizes.



Downhill walking as an option (-10%)

This treadmill can be executed with 10% negative elevation. This allows for downhill walking which is extremely useful for rehabilitation of certain injuries.



Optional integrated BPM and SpO2

The treadmill can be extended with a stable and reliable blood pressure module and SpO2 measurement.



LCRM compatible

This treadmill can be connected to the Lode Cardiac Rehab manager offering a complete cardiac rehab setting

17
cm

Extreme low step-up height

To allow people to safely and comfortably step up the treadmill it is important to have a very low step-up height

0.1
km/h

Small adjustment steps

The speed of the treadmill can be adjusted in the smallest steps of only 0.1 km/h!



Low cover plate

The motor compartment of the treadmill is designed in such a way that the cover is only marginally higher than the belt surface. This allows for low camera positions in case the treadmill is used for gait analysis purposes. Also a therapist has the best possible view on feet and lower extremities.

Valiant 2 rehab XXL

Robust and reliable treadmill for rehabilitation



Interconnectivity between Lode products

Connecting Lode products has never been easier! Lode rehab and sports products have a standard Network card:

- To be able to connect the first product to the PC with L(C)RM a Lode proprietary network to PC cable is needed (#930930). This cable is standard included with Lode Rehab Software.
- From the second product onward products can be connected to the previous one, creating a bus network configuration;
- The last product always needs a termination plug to avoid interference and loss of data. Therefore all products with such a network card come with a termination plug.

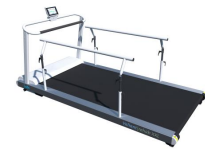


Benefits



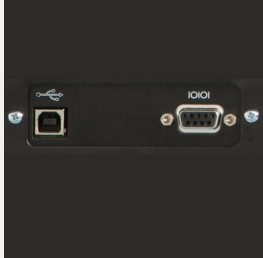





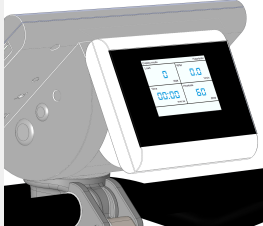
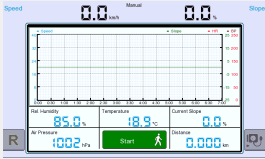




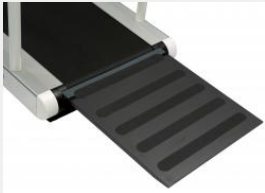
- Lossless data connection
- High bandwidth
- No interference of COM ports
- Daisy chain connection
- Full access of all data in the product to LCRM

Valiant 2 rehab XXL

Robust and reliable treadmill for rehabilitation



Valiant 2 rehab XXL can a.o be extended with the following options:

<p>USB to Serial converter</p> <p>Easy connection</p>  <p>Partnumber: 226012</p>	<p>Emergency Stop Button</p> <p>Ultimate safety</p>  <p>Partnumber: 945804</p>	<p>Communication Module</p> <p>Connect to ECG and spirometers</p>  <p>Partnumber: 945850</p>	<p>Extension for emergency lanyard</p> <p>Original Accessory</p>  <p>Partnumber: 945931</p>	<p>Blood Pressure Measurement with ECG trigger for treadmills with ECG trigger</p>  <p>Partnumber: 945824</p>
<p>SpO2 for control unit with touch panel (extra long cable)</p> <p>Oxygen saturation</p>  <p>Partnumber: 945822</p>	<p>Handrails, Side - Adj. for extended Valiant 2</p> <p>Safe and flexible handrails</p>  <p>Partnumber: 938821</p>	<p>Table foot for Control Unit with Touchscreen</p> <p>Control from any desired place</p>  <p>Partnumber: 945818</p>	<p>Colour Display 3.5" - 2nd screen</p> <p>Multifunctionality</p>  <p>Partnumber: 945819</p>	<p>Ambient sensor pack</p> <p>Check environmental conditions during test</p>  <p>Partnumber: 945826</p>
<p>Network module connection cable</p> <p>Easy connection</p>  <p>Partnumber: 930930</p>	<p>SpO2 for control unit with touch panel - ordered afterwards</p> <p>Ordered afterwards</p>  <p>Partnumber: P945822</p>	<p>Add program function to 7" touch screen for Treadmill</p> <p>Programmable</p>  <p>Partnumber: U945815</p>	<p>Remote Control for treadmill</p> <p>Ease of use</p>  <p>Partnumber: 945831</p>	<p>Entrance plate for Valiant rehab XXL</p> <p>Even easier entrance to the treadmill</p>  <p>Partnumber: 938849</p>

Valiant 2 rehab XXL

Robust and reliable treadmill for rehabilitation



Specifications

Workload

Maximum speed	20 km/h	12.4 mph
Minimum operational speed	0.1 km/h	0.1 mph
Positive elevation	15 %	
Elevation adjustment steps	0.5 %	

Accuracy

Speed accuracy	5 %	
Accuracy inclination	0,5 %	

Comfort

Allowed user weight	225 kg	496 lbs
---------------------	--------	---------

User Interface

English user interface	✓
Norwegian user interface	✓
Czech user interface	✓
Danish user interface	✓
Dutch user interface	✓
Finnish user interface	✓
French user interface	✓
German user interface	✓
Italian user interface	✓
Japanese user interface	✓
Korean user interface	✓
Polish user interface	✓
Portugese user interface	✓
Russian user interface	✓
Spanish user interface	✓
Turkish user interface	✓
Ukrainian user interface	✓

Connectivity

Optional USB connector	✓
Optional RS232 connector	✓

Dimensions

Screen resolution	800 x 400 pixels	
Walking surface length	250 cm	98.4 inch
Walking surface width	120 cm	47.2 inch
Step up height	17 cm	6.7 inch
Product length (cm)	312 cm	122.8 inch
Product width (cm)	148 cm	58.3 inch
Product height	130 cm	51.2 inch

Order info

Partnumber: 938915

*Specifications are subject to change without notice.

Power requirements

Power cord length	250 cm	98.4 inch
Power cord IEC 60320 C19 with CEE 7/7 plug	✓	
Power cord NEMA	✗	
Maximum rated current	12 A	
115 V AC 50/60 Hz (2 phases)	✓	
230 V AC 50/60 Hz	✓	
Maximum motor power	3.0 kW	

Standards & Safety

IEC 60601-1:2005	✓
ISO 13485:2016 compliant	✓
ISO 9001:2015 compliant	✓
Standard emergency lanyard	✓